

Auto Fuser K Series

Disposable elastomeric infusion pump



Auto Fuser K

Elastomeric infusion pump is designed to provide clinicians with the option of delivering pre-determined volume of medication at continuous and accurate flow rate for patients at the hospitals, ambulatory settings or homecare.

• Elastomeric infusion pump

Easy to check and protect the remaining medication with transparent and protective housing

• Fill port

Designed to prevent reflux or leakage while filling




• Clip

Can be attached to clothes, pillow or shower curtain to prevent pumps from falling




Luer lock Syringe Connecting

The luer lock type filling port, designed to prevent leakage in the filling process and the flat surface of the pump's bottom allows two different ways of filling; as shown in the picture.

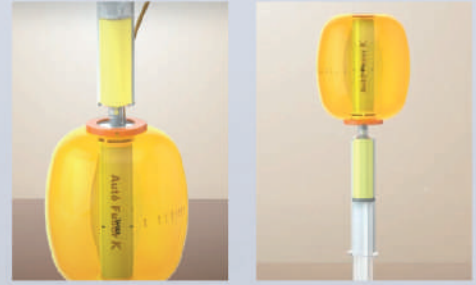

• Air filter

Filtering air bubble and foreign particles

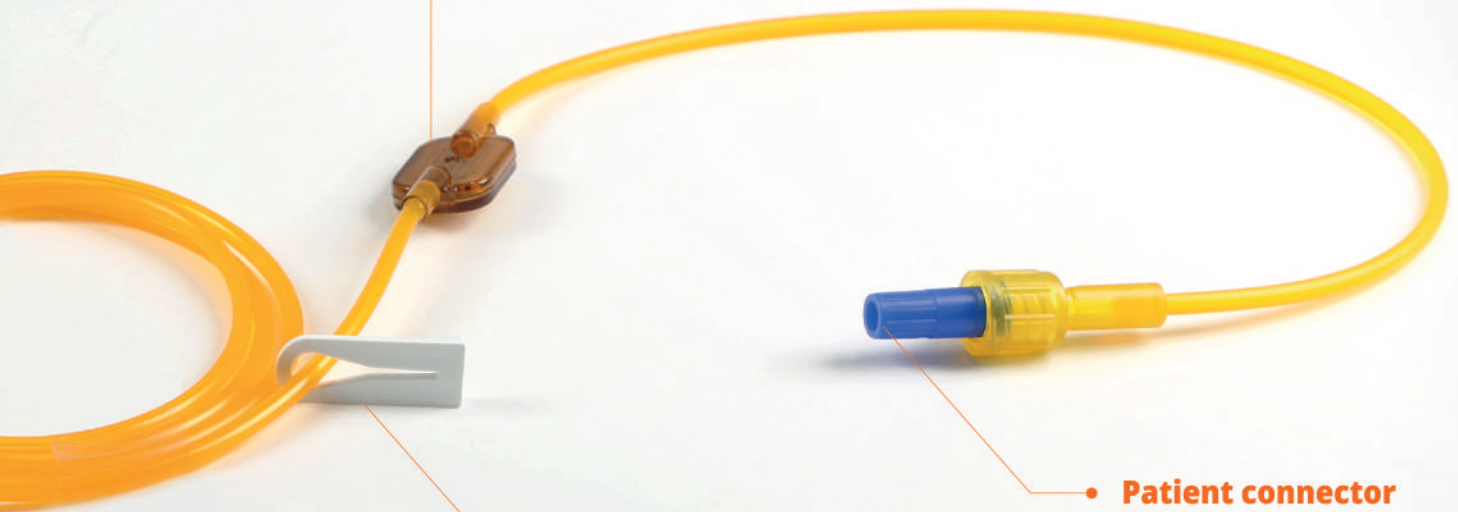

• Patient connector

Convenient priming due to purge filter

• Clamp

Preventing kinking and any damages on tubing due to its meticulous and blunt design

Auto Selector K

*AutoSelector K provides selectable flow rate
with dial type flow regulator
It can combine with Bolus button for PCA application*



Selector

- Dial can be removed so the patient cannot temper with the flow rate, only the medical staff can handle setting the flow rate
- Diverse flow rate selectable (2 types of Selector)

Selective Flow Rate Infuser

Selector can be adjusted the flow rate in 8 discrete steps from 0.5 ml/hr – 7 ml/hr or 1 ml/hr – 14 ml/hr according to your patient's needs. Multiple choices in flow rates and PCA options allow you to customize your patient's pain protocol.



A Type



B Type



Bolus


- Easy and convenient use due to its pen type design
- Bolus dose can be observed through window in real time



Advantages

Compact Size

Minimizing the dead space and unnecessary space between silicone balloon and hard case, AutoFuser K is compact, which allows better portability

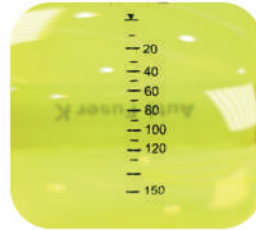

Hardcase Design

Hard case protects inner elastomeric membrane against physical impact. The housing is completely transparent, which allows clinicians and patients to check remaining medication and its proper operation.

UV Protection

Medication inside the pump can be protected from UV light with UV protection material.



Easy Filling

One hand filling is available, which provides medical/nursing staff with better efficiency

Patient Connector

Convenient priming process due to purge filter at patients end



Latex free & DEHP free

AutoFuser K is latex free and DEHP free, prevents allergic symptoms for patient safety

Indications



Post-operative pain management



Cancer pain management



Chemotherapy



Thalassemia




Antibiotics



Hydration

Product Portfolio

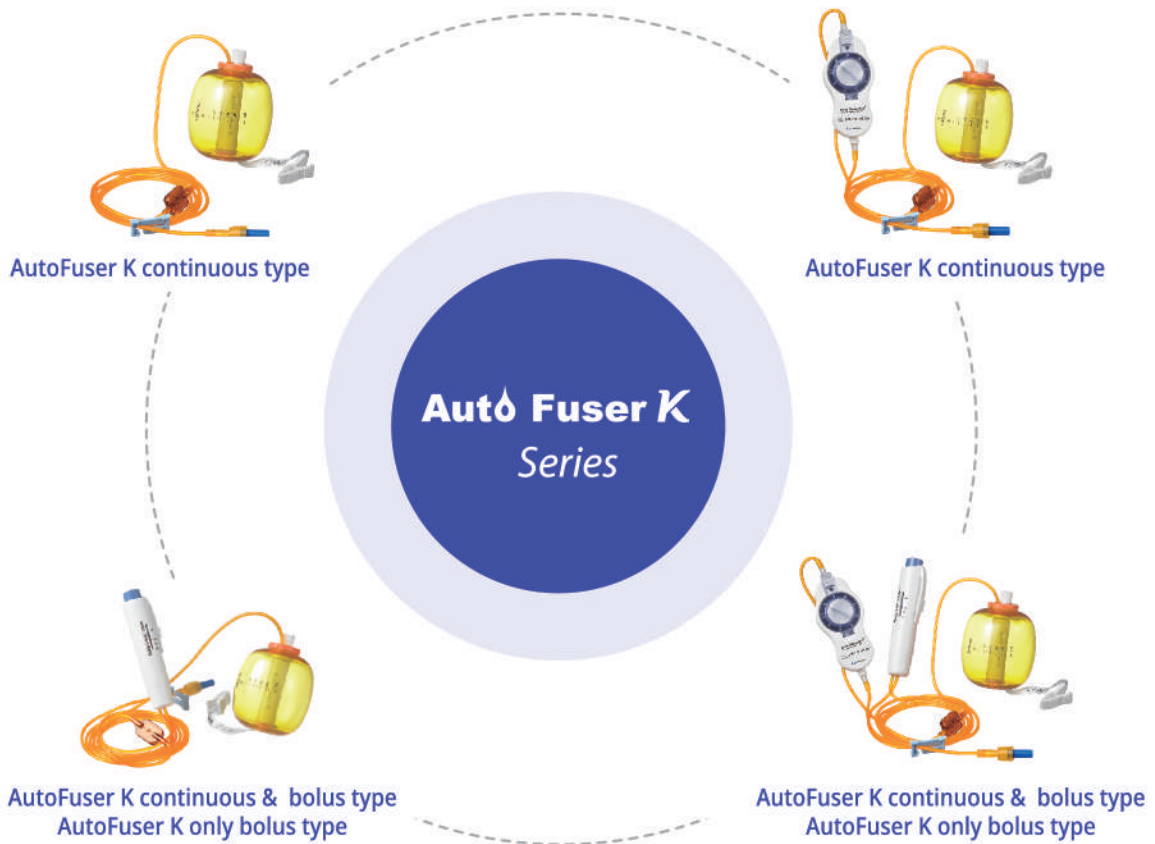

Specification



Volume size	60, 100, 150, 250, 400, 500mL
Flow rate	0.5 ~ 500mL/hr
Selector rate	0.5 to 7 mL/hr (A-type) 1 to 14 mL/hr (B-type)
Bolus Volume	0.5, 1.0, 2.0, 3.0, 5.0mL
Lock out Time	0, 15, 30, 60min

Auto Fuser *K* Series

ACEMEDICAL

33, Naeyu-gil, 124beon-gil, Deogyang-gu, Goyang-si, Gyeonggi-do, Korea
Tel : +82-2-920-0001 Fax : +82-2-922-9319 <http://www.ace-medical.com>